RELIABLE BEAUTY SOLUTIONS

PiiGii

PG-BD05

Lumi Tight Beauty Device

User Manual

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2-Year Replacement Service

Members can upgrade to a 2-year replacement service by registering via WeChat or WhatsApp

— SCAN NOW TO UPGRADE —

Thank you for choosing PiiGii.

To ensure safe and proper use, please read this user manual carefully before using the product and follow all instructions. Keep this manual for future reference. If you give this product to someone else, please make sure to pass this manual to the new user as well.

Table of Contents

Safety Precautions	34
Before First Use	39
Product Overview and Accessories	40
Operating Instructions	41
How to Use	43
Product Care and Maintenance	45
Product Specifications	46
Troubleshooting	46
Warranty Information	47
Disposal Statement	48
Disclaimer	48

SAFETY PRECAUTIONS

- **WARNING:** Failure to follow these instructions may result in personal injury or property damage.
- Please read this Safety Precautions section carefully before using the product and follow all instructions.
- Use the product only as directed and in accordance with the operating procedures in this manual.
- If anything feels wrong or unusual during use, stop using the device immediately.
- Using the device while it is damaged or not working properly may cause personal injury or property damage.
- To help prevent personal injury or property damage to you or others, this manual uses the following terms to indicate different levels of risk:

⚠ DANGER: Incorrect use will result in death or serious injury.

⚠ WARNING: Incorrect use could result in death or serious injury.

⚠ CAUTION: Incorrect use may result in injury or property damage.

🚫 PROHIBITED: Indicates actions that must not be performed.

⚠ DANGER

Do not use this device while wearing or using any other electronic device. Do not use this device with any medical electronic equipment, as this may cause the equipment to malfunction and result in serious injury or health problems.

This includes, but is not limited to:

- Pacemakers or other implanted medical devices.
- Life support equipment such as heart-lung machines.
- Wearable medical devices such as ECG monitors or IV equipment.

⚠ WARNING

Never use this device if you:

- Are currently receiving medical treatment, have a fever, or have an infectious or acute illness.
- Have a serious internal organ disease.
- Have tuberculosis.
- Have heart disease, abnormal blood pressure, or cancer.
- Have hemophilia, any bleeding disorder, diabetes, or facial neuralgia.
- Have a weakened immune system due to HIV or other immune disorders, or are using immunosuppressant medicines.

- Are currently receiving dental treatment, or have had dental restoration work in the treatment area.
- Have abnormal sensation or pain, or cannot properly feel heat or cold.
- Tend to develop keloid scars, or have facial hyperpigmentation, photosensitivity, or other pigment disorders.
- Take steroid medicines for a long time, or have dilated capillaries due to liver problems.
- Have allergies, hereditary or atopic dermatitis, other skin diseases, very sensitive skin, or very fragile skin.
- Have a history of epilepsy, or feel unwell when you see flashing lights.
- Have after effects in the head or eyes due to a traffic accident or other injury.
- Have had surgery on your face or eyes, including laser eye surgery such as LASIK.
- Are an infant, a young child, or cannot clearly express your own feelings or discomfort.
- Are under the age of 16 years.
- Are highly short sighted (more than -6.0 diopters).
- Are pregnant, on your period, or breastfeeding.
- Are drunk or hungover.
- Are taking sleeping pills.
- Have recently used medicines that increase sensitivity to light, such as tetracycline.
- Are using immunosuppressants, or are taking medicines such as aspirin or isotretinoin.
- Have any bleeding disorder.
- Are using topical products containing isotretinoin, alpha hydroxy acids (AHA), beta hydroxy acids (BHA), or azelaic acid.
- Have used whitening or brightening skincare products in the past week, or have used cosmetics that contain vitamin A.
- Have skin that has been damaged by sun or other light exposure.
- Are extremely tired or feeling unwell.
- Have had pigmentation, dermatitis, or other skin problems caused by synthetic fibers.
- Have had dermatitis or other skin problems caused by cosmetics.
- Have pigmentation caused by physical irritation.
- Have inflamed, red, or swollen skin.
- Have skin with pus or acne.
- Have recently received cosmetic or aesthetic treatments at a beauty salon or clinic, such as peeling treatments.
- Have had skin problems in the past after using beauty devices.

※ If you have any other medical condition, are under medical care for any reason, or are not sure whether this device is safe for you, do not use it until you have consulted your doctor or healthcare professional.

PROHIBITED AREAS AND SKIN CONDITIONS

Never use this device on the following areas or skin conditions:

- Around the heart or directly over major blood vessels, such as large veins and arteries.
- The eyeball, upper eyelids, eyebrows, scalp, or ears.
- The lips, inside the mouth, the thyroid area, or around the Adam's apple.
- The nipples, areola, genital area, anus, or any other mucous membranes.
- Areas with conditions such as melasma, abnormal dark spots, moles, or viral warts.
- Areas that feel tingling, red, inflamed, itchy, hot, or otherwise abnormal.
- Areas with wounds, eczema, psoriasis, herpes, pustular acne, or severe pimples.
- Areas that have had surgery and are still healing.
- Areas that have recently been treated with hair removal, waxing, or shaving.

PROHIBITION

- Do not use this device together with any other device at the same time.
- Do not use the device continuously for a long time, as this may cause burns.
- Do not use the device on any skin areas where the treatment head cannot fully touch the skin.
- Do not use the device for a long time on areas with thin skin or directly over bones.
- Always move the device over the skin. Do not hold it in one spot for a long time.
- Do not repeatedly treat the same area. This may damage your skin or adversely affect your health.
- Do not use the device if it is damaged or not working properly.
- Do not use the device for any purpose other than cosmetic skincare.
- This device is for beauty care only. Do not use it for any treatment purposes or for any other use not described in this manual.
- Do not use the device on animals or pets.
- Do not burn or incinerate the device, and do not use or store it in high temperature environments, as this may cause overheating, fire, or malfunction.
- Do not use or store this device in bathrooms, shower rooms, bathtubs, washbasins, on washstands, or in any other damp, humid, or splash-prone places.
- Do not immerse this device in water or allow water to splash or drip onto it, as this may cause electric shock, short circuit, fire, or malfunction.
- Do not touch the device with wet hands, as this may cause electric shock, short circuit, fire, or malfunction.
- Do not use the device near flammable materials such as alcohol, thinners, solvents, nail polish remover, or sprays, as this may cause overheating, fire, or malfunction.

- Do not use the device in very dusty places or in direct sunlight.
- Do not use the device in places with strong electromagnetic waves or strong magnetic fields.
- Do not put the device into a microwave oven or any high pressure container or vessel.
- Do not use the device in places with sudden changes in temperature, such as right next to air conditioners or heaters.
- If you notice a strange smell or abnormal heat from the device, stop using it immediately, as continued use may cause damage or fire.
- Do not place heavy objects on the device, and do not drop or hit the device during use, as this may cause malfunction or injury.
- Do not place the device on unstable, sloping, vibrating, or impact-prone surfaces.
- Do not allow people who cannot clearly express their own intentions to use this device, and avoid use by people with disabilities. This may cause accidents. If necessary, consult your doctor or healthcare professional before use.
- People under 16 years of age must not use this device.
- Children must not use, clean, or service this device without supervision.
- People with reduced physical, sensory, or mental abilities, or those who lack experience and knowledge, must not use this device unless they are supervised or instructed.
- Do not use this device on broken skin.
- Use only disposable and consumable parts recommended by the manufacturer.
- Do not modify or disable any of the built in safety features of this device.

CAUTION

- Do not press or hold the device on the same spot, keep it moving to avoid skin injury.
- Do not allow the device to come into contact with or strike other objects, as this may cause overheating, fire, malfunction, or damage.
- Before use, check the treatment head for cracks, damage, or deformation. If you notice any of these problems, do not use the device, as continued use may cause skin injury or malfunction.
- Do not use the device in extremely hot or extremely cold environments.
- The recommended operating temperature is about 15°C to 30°C.
- Do not leave the device on the floor or on furniture while it is switched on.
- Turn off the power immediately after use to help prevent damage to furniture or injury.
- While the power is on, do not use the device in contact with other electrical products such as electric blankets, as this may cause malfunction or injury.
- Do not use this device in contact with or very close to metal objects or other electrical products, such as mobile phones or televisions.
- Do not use this device together with exfoliating skincare products, as this may cause skin problems.

- Do not look directly at the light for a long time, as this may harm your eyes.
- Regularly check the device and gently clean the glass area to remove dirt, to help prevent malfunction or damage.
- Do not use the device if the air vents are blocked, as this may cause malfunction or fire.
- Do not shine the light on clothing, hair, or other objects, as this may cause scorching, discoloration, or damage.
- If you have recently received cosmetic or aesthetic treatments, consult your doctor or healthcare professional before using this device.
- If your skin feels slightly warm or looks slightly red after use, this is normal and can be relieved with moisturizing and hydrating skincare.
- If the pain is strong, select a lower EMS level or stop using the device.
- When treating areas close to the eyes, wear suitable protective goggles and keep your eyes closed.
- Never leave the device switched on and unattended. After use, switch it off, unplug it, and store it properly.
- For 48 hours after treatment, avoid exposing the treated area to direct ultraviolet light.
- Do not leave the device in direct sunlight or under ultraviolet light for a long time, as this may seriously damage it.
- Do not attempt to disassemble, repair, or modify the device, power adapter, or power cord yourself.
- Do not bend, pull, or twist the power cord.
- Do not use or store the device in water, in bathrooms, on washstands, or in other very humid places.
- Keep the device, power adapter, and power cord away from water.
- Do not use power adapters or power cords that are not supplied or approved by the manufacturer.
- Because plug shapes and voltages may differ by region, always check that they are suitable before use.
- If the power supply is unstable or outside the range of 100 to 240V AC, use a voltage stabilizer or transformer before operating the device, otherwise the performance of the device may be affected.
- Do not connect this product in parallel or in series with other electrical appliances.
- Using multiple devices in parallel beyond the rated load may cause fire.
- If a power outage occurs during use, unplug the power cord from the outlet.
- Use only the power adapter supplied with this product. Using an unsuitable adapter may damage the device.

PROHIBITED DISASSEMBLY

※ Do not disassemble, modify, or repair this device yourself.

If the device needs inspection or repair, please contact our customer service.

CAUTIONS FOR USE


- This device may become warm during use.
- Before each use, clean the treatment head.
- Make sure the treatment area and the surrounding skin are free of jewellery, glasses, contact lenses, hair, and any other objects or devices before using the device.
- When using the device for the first time, or if you are not yet used to it, always start with the lowest EMS level.
- Do not press the treatment head hard on your skin, and do not keep it on the same spot for more than three seconds.
- Do not use the same mode more than once a day, as overuse may cause skin problems.
- When treating the area around the eyes, wear suitable protective goggles and keep your eyes closed.
- During use, apply water based skincare products as needed to keep your skin well moisturised.
- If you feel strong stinging on your skin, select a lower EMS level or move to another area. Apply gel or other water based skincare products before continuing.
- If you feel pain or any unusual discomfort in your skin or body, stop using the device immediately and consult your doctor or healthcare professional.
- Do not treat the skin too often. Use the device according to the condition of your skin and follow the recommended frequency of use.
- For facial care, remove makeup, wash your face, and then apply water based skincare products before using the device.
- If your skin develops a rash, redness, itching, or any other abnormal reaction during or after use, stop using the device and consult your doctor or healthcare professional.
- After use, do not touch the glass area of the device immediately, as it may still be hot.

BEFORE FIRST USE

- After opening the box, take out the main device and all accessories. Keep the packaging and this manual for future reference.
- Check whether the device and all accessories are in good condition. If you find any damage or the device is not working properly, please contact the distributor or our customer service.
- If the light intensity of the device is unstable or abnormal, switch off the device immediately and contact the distributor or our customer service. Before use, clean the skin and make sure it is dry. Always use this device together with a gel.

PRODUCT OVERVIEW AND ACCESSORIES

Names of Parts

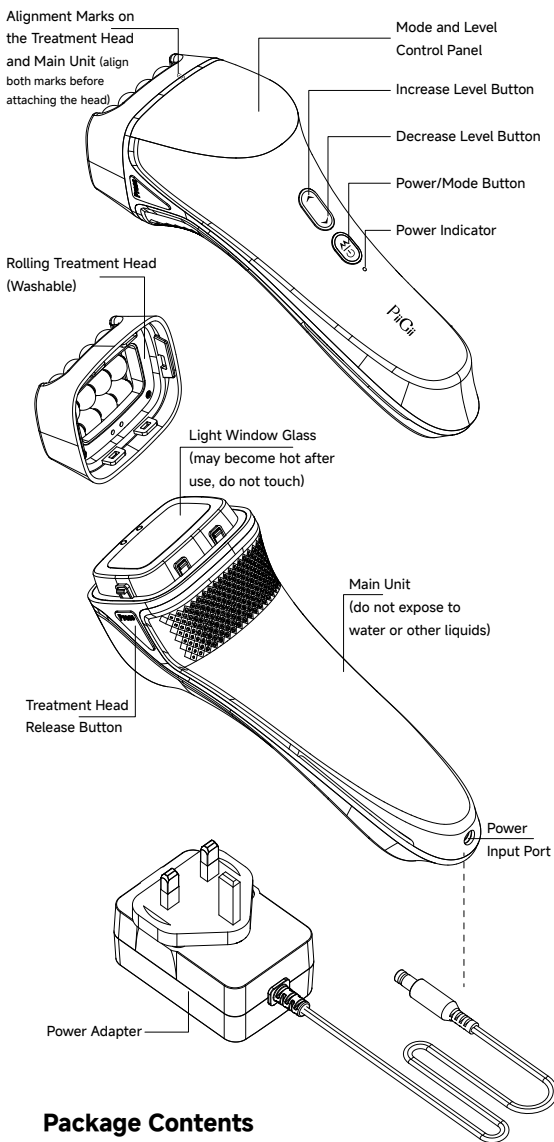Package Contents

Main Unit ×1

Rolling Treatment Head ×1


Power Adapter ×1

User Manual ×1

Quick Start Guide ×1

OPERATING INSTRUCTIONS

- 1. Connect the power:** Connect the power adapter to the device, then plug the adapter into a wall socket. The power indicator will turn red.
- 2. Selecting a Mode:** Press and hold the "Power/Mode Button" for 2 seconds to switch on the device. The power indicator will change from red to white and the default Face Mode will turn on. Short press the "Power/Mode Button" to switch between the 5 modes: Face, Waist, Arms, Legs and Bust. The indicator for the selected mode will turn on.

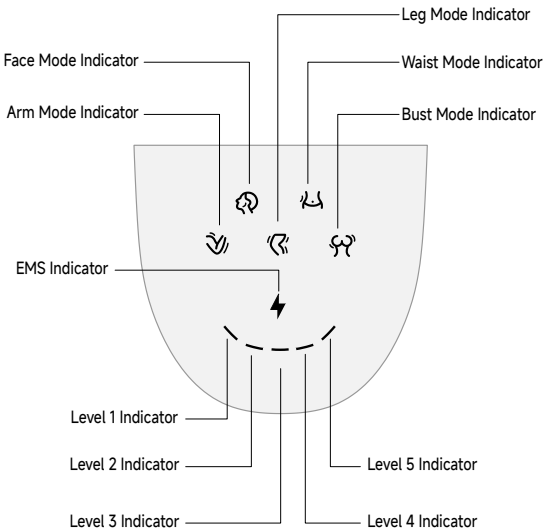- 3. Set the EMS level:** The EMS function is off by default. After selecting a mode, press the "Increase Level Button" to turn on EMS. The level will go up from 1 to 5 and the EMS indicator and the selected level indicator will light up. Press the "Decrease Level Button" to lower the EMS level step by step from 5 back to 1. When EMS is at level 1, press the "Decrease Level Button" again to switch EMS off.

Note: Higher EMS levels mean stronger intensity. In Bust mode, the EMS function is not available.

- 4. Start the treatment:** After a mode is selected, the red light will turn on and the device will start working. If the device is kept still for about 3 seconds, the red light will turn off automatically but EMS output will continue. When you move the device again, the red light will turn back on. When the red light is on, the built in cooling fan will also start. After the red light turns off, the fan will continue running for a short time and then switch off automatically.

- 5. Treatment time:** The preset treatment time for each mode is 10 minutes. After 10 minutes of use, the device will automatically end the current mode, beep once, and enter standby mode. At this time, the current mode indicator will flash, and the red light and EMS functions will stop immediately, while the fan will turn off automatically after a short delay. If no mode is selected within 3 minutes in standby mode, the device will power off automatically.


- 6. Switch modes:** Switch modes: During treatment or after a treatment ends, you can short press the "Power/Mode Button" to switch to the next mode. The 10-minute timer will restart for the new mode.

- 7. End the treatment:** During treatment or after a treatment ends, press and hold the "Power/Mode Button" for 2 seconds to switch off the device. All functions will stop immediately, the power indicator will change from white to red, and the fan will turn off automatically after a short delay.

- 8. After switching off:** Wait until the device and treatment head are completely cool. Then unplug the power adapter from the wall socket and remove the treatment head for cleaning.
Caution: During and after use, the glass area of the Light Window Glass can become very hot. Touching it may cause burns. Do not touch this area.

HOW TO USE


01 Install and remove the treatment head

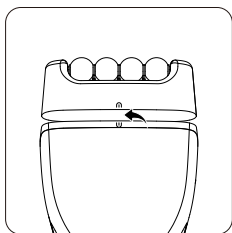INSTALL

Align the marks on the treatment head and the main unit. Gently press the treatment head down until it clicks into place.

Note: The treatment head is already installed when you receive the device.

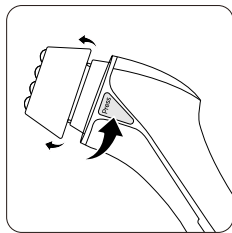REMOVAL

Press the release buttons on both sides and pull the treatment head upwards to detach.

02. Treatment steps

1 Remove makeup and cleanse

Before use, remove makeup and wash your face. Take off any accessories or items on the face, ears, neck or other treatment areas, including contact lenses. Not removing these may cause injury or device malfunction.


2 Apply gel

For best results, use the device together with the included gel. Apply an even layer of gel to the face. If your skin feels dry or the treatment head does not glide smoothly during use, apply more gel as needed so the head can move smoothly over the skin.

③ Switch on and use

After connecting the power adapter, press and hold the "Power/Mode Button" for 2 seconds to switch on the device. Select the desired mode and set the EMS level. The red light will turn on automatically. Make sure the treatment head is in full contact with the skin and then start moving it. **Note:** For first time use, always start from EMS level 1.

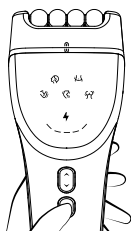④ Treatment method


Gently press the treatment head and roll it back and forth over the treatment area to help your skincare products absorb better. The preset treatment time is 10 minutes. After 10 minutes, the device will stop working, beep once, and enter standby mode.

⑤ Clean and care for the skin


After the treatment, rinse your face with clean water. Then follow your usual skincare routine and apply suitable products for hydration and moisture.


03 Using on other body areas

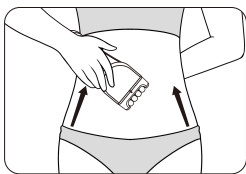
Waist and abdomen

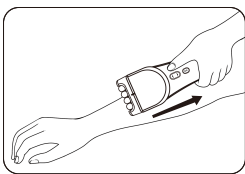Arms

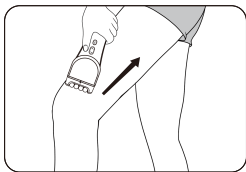Legs


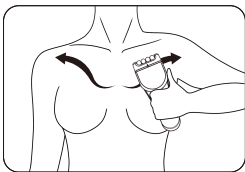
Bust

Important notes:

- Modes and intensity levels are for intended body areas only. Use only as directed. Improper use may result in injury.
- Use the device directly on the skin. Do not use it over clothing.
- Always use the device with water-based skincare products.
- Do not use oily or scrub type products, as they may cause device damage or skin problems.
- The above methods are for reference only. You may adjust according to your personal needs.

04 Cleaning the device

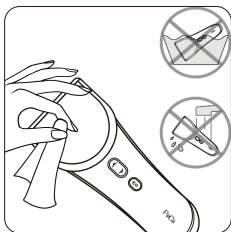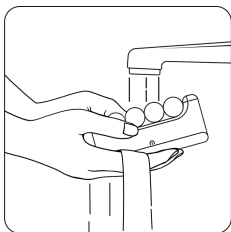※ **Important:** After the treatment, do not touch the glass area on the main unit. During and after use, the glass area can become very hot and may cause burns. Always wait until the device and treatment head are completely cool before removing the head.

- Before cleaning the device, switch it off, unplug the power adapter and allow the device to cool down.
- Use soft tissues or purified water wipes to wipe the glass area and the surface of the device. Clean the glass area after each use.
- If tissues or purified water wipes are not enough to clean the glass area, you may use wipes with 75 percent medical alcohol. Do not use steel wool, abrasive tools or corrosive cleaners.
- The main unit is not waterproof. Do not wash it or put it in water. Only the treatment head is washable.
- This device is intended for personal use. Do not share it with others. If sharing is necessary, always clean and disinfect the glass area before each use to reduce the risk of cross infection.

PRODUCT CARE AND MAINTENANCE

- Keep the device dry and away from water to avoid the risk of electric shock.
 - To avoid damage or loss of the device and its accessories, store them in the box after each use.
 - If you will not use the device for a long time, store it properly in a dry, cool place away from direct sunlight. Avoid contact with hard objects to prevent damage.
 - Do not store the device or its accessories in a damp or humid environment, as this may cause malfunction.
 - Do not clean the device while it is plugged in.
 - Before each use, check the device and the power cord for any signs of damage. If you find any damage or cracks, stop using the device immediately and contact our customer service.
 - Always unplug the power adapter from the wall socket when the device is not in use.
- ※ If the device has not been used for six months or longer, check it carefully and make sure it is working normally before use. Do not use the device if it is damaged or malfunctioning, as this may result in injury.

PRODUCT SPECIFICATIONS

Product Name	Lumi Tight Beauty Device
Model No.	PG-BD05
Power Supply	AC Adapter (RY36A120300BS)
Cable length	1.5 m
Adapter Input	AC 100-240V, 50/60Hz 1.0A
Adapter Output	DC 12.0V \equiv 3.0A
Standards	GB 4706.85-2008, GB 4706.15-2021
	GB 4706.1-2024, GB/T 36419-2018
Dimensions	202 × 70 × 82mm
Weight	Approx. 280g

TROUBLESHOOTING

If the following steps do not solve the problem, please contact our customer service team.

Problem	Possible Cause	Solution
Unable to power on	Power not connected	Check if power is connected
No light output	The device or treatment head is too hot	Wait until the device has cooled down, then try again
Power indicator does not turn on	Poor power connection	Check the power connection
Skin becomes red or sensitive after use	Temporary redness or sensitivity, which is usually normal	Stop using the device and wait until the skin calms down
Unsatisfactory results after use	The wrong mode was selected for the treatment area	Select the mode that matches the treatment area
	Sun protection after treatment was not enough	Use proper sun protection after each treatment
Device becomes hot during use	The device stayed too long on the same area	Move the treatment head more quickly and avoid staying on one spot
	The device has been used continuously for a long time	Switch off the device and let it cool down before continuing

WARRANTY INFORMATION

Standard Warranty: 2 years

Register via WeChat or WhatsApp to upgrade to a 2-year replacement service.

The free warranty service does not cover:

- Any damage caused by not following the user manual.
- Damage or loss caused by human error, misuse, abuse, negligence, accidental damage, dropping, rust, natural disasters, or any other causes beyond our control.
- Problems caused by incorrect installation of parts or accessories.
- Use of the product for commercial purposes or for any purpose that is different from its original intended use.
- Any unauthorised repair or modification.
- Normal wear and tear (for example, scratches, paint chips and other cosmetic issues).
- Any signs of moisture, water or liquid damage.

Our company is not responsible for any indirect or consequential loss caused by the use of this product. This includes but is not limited to loss of business income, business interruption, the cost of renting or buying replacement equipment, and any claims made by third parties. **For full warranty terms and conditions, please visit our official website.**

Hong Kong Customer Service Center

Regal Silver Limited.

Tel Number: +852-2186 6089

Address: Unit 1307, Vanta Industrial Centre, 21-33 Tai Lin Pai Road,
Kwai Chung, N.T., Hong Kong

Email: cs@piigii.com

WhatsApp: +852-9806 1578

Website: <http://www.piigii.com/zh-hk>

Shenzhen Customer Service Center

Shenzhen Tian Hao Regal Silver Technology Limited.

Tel number: +86-1337 866 5715

Address: No. 181, Jiangshi Road, Shiwei Community, Matian Street,
Guangming District, Shenzhen City, Guangdong, China

Website: <http://www.piigii.cn/>

For inquiries regarding product maintenance, please contact the customer service center.

DISPOSAL STATEMENT

When disposing of this device, do not throw it away with household waste. Please follow local laws and regulations for waste and electronic products. Proper disposal helps protect the environment.

DISCLAIMER

Regal Silver Limited has made every effort to ensure that the information in this manual is current and accurate. However, the company is not responsible for any errors or omissions, and the contents are subject to change without prior notice.

Regal Silver Limited is not liable for any loss or damage arising from misuse of the product, use that does not follow this manual, or use in any way other than the methods described in this manual, on the product packaging or in promotional materials. When using this device, users are responsible for reading and following all safety information.

This manual may not be reproduced, stored or translated in any form or by any means without the written authorization of Regal Silver Limited.

Regal Silver Limited does not guarantee the results of using this device. Results may vary depending on the user's individual skin condition. It is strongly recommended that all users consult your doctor or healthcare professional before using any skincare or beauty device.

These statements have not been evaluated by the Food and Drug Administration. This product is not intended to diagnose, treat, cure or prevent any disease. If you have any questions or health concerns, please consult your doctor or healthcare professional before use.

PiiGii

www.piiyii.cn

www.piiyii.com/zh-hk