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
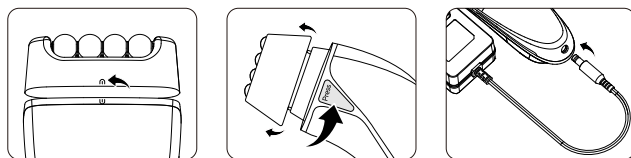
## PRODUCT INSTALLATION
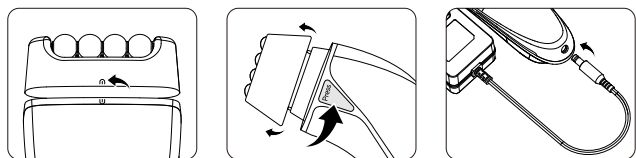
### INSTALL

Align the marks on the treatment head and the main unit. Gently press the treatment head down until it clicks into place. Note: The treatment head is already installed when you receive the device.

### REMOVAL

Press the release buttons on both sides and pull the treatment head upwards to detach.

### POWER CONNECTION

Connect the output connector of the power adapter to the main unit, and insert the adapter plug into the power socket. When the power indicator on the main unit lights up in red, the device is in standby mode.

## TREATMENT STEPS

### 1 Remove makeup and cleanse

Before use, remove makeup and wash your face. Take off any accessories or items on the face, ears, neck or other treatment areas, including contact lenses. Not removing these may cause injury or device malfunction.


### 2 Apply gel

For best results, use the device together with the included gel. Apply an even layer of gel to the face. If your skin feels dry or the treatment head does not glide smoothly during use, apply more gel as needed so the head can move smoothly over the skin.

### 3 Switch on and use

After connecting the power adapter, press and hold the "Power/Mode Button" for 2 seconds to switch on the device. Select the desired mode and set the EMS level. The red light will turn on automatically. Make sure the treatment head is in full contact with the skin and then start moving it. Note: For first time use, always start from EMS level 1. Modes and intensity levels are for intended body areas only. Use only as directed. Improper use may result in injury.
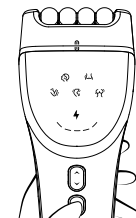
### 4 Treatment method

Gently press the treatment head and roll it back and forth over the treatment area to help your skincare products absorb better. The preset treatment time is 10 minutes. After 10 minutes, the device will stop working, beep once, and enter standby mode.


### 5 Clean and care for the skin

After the treatment, rinse your face with clean water. Then follow your usual skincare routine and apply suitable products for hydration and moisture.

## CLEANING

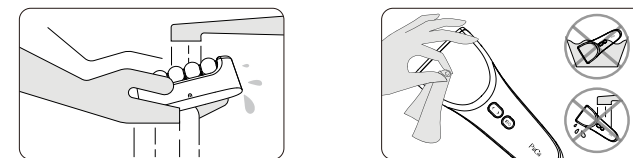**⚠ Important:** After the treatment, do not touch the glass area on the main unit. During and after use, the glass area can become very hot and may cause burns. Always wait until the device and treatment head are completely cool before removing the head.

- Before cleaning the device, switch it off, unplug the power adapter and allow the device to cool down.
- Use soft tissues or purified water wipes to wipe the glass area and the surface of the device. Clean the glass area after each use.
- If tissues or purified water wipes are not enough to clean the glass area, you may use wipes with 75 percent medical alcohol. Do not use steel wool, abrasive tools or corrosive cleaners.
- The main unit is not waterproof. Do not wash it or put it in water. Only the treatment head is washable.

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
## MEMBER REWARDS

JOIN AS A MEMBER BY SCANNING THE QR CODE

—— ENJOY THESE BENEFITS ——

- ① TWO YEARS FREE REPLACEMENT
- ② PERSONAL SKINCARE ADVISOR
- ③ SPECIAL MEMBER REWARDS
- ④ DEVICE USAGE GUIDANCE

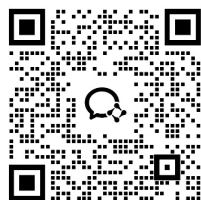\* REGISTER WITHIN 1 MONTH OF PURCHASE TO UNLOCK THE 'TWO YEARS FREE REPLACEMENT' SERVICE

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